

Mattresses

Options available may include Foam, Latex or specifically designed Inner spring. Conventional inner spring mattresses are not suitable for electrically adjustable beds.

See AT Australia guidelines for 'Choosing a Bed and Mattress.'

Other Options

To increase the height of the bed, use 'bed blocks'.

To elevate the feet, use a wedge-shaped cushion.


To get in and out of bed, try bed rails.

To sit up supported in bed, options include:

- an angled backrest frame, which sits on top of the mattress, to support the pillow
- an electrically adjustable backrest frame which sits under the mattress
- a 'bed lounge chair' foam support

Assistive Technology Australia
Shop 4019, Level 4
Westpoint Shopping Centre
17 Patrick Street
Blacktown NSW 2148

(previously known as Independent Living Centre NSW)

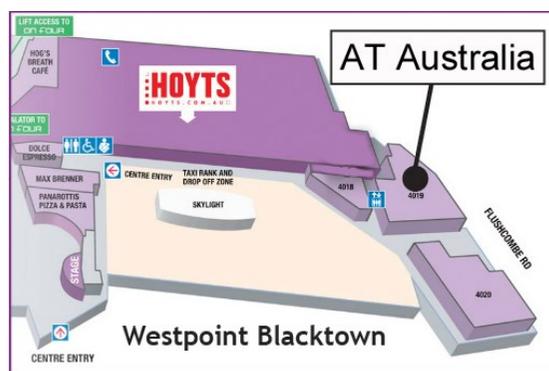

Info-Line: 1300 452 679
Reception: (02) 9912 5800
Fax: (02) 8814 9656
Email: help@at-aust.org
Web: www.at-aust.org

Supported by



**Family &
Community Services**
Ageing, Disability & Home Care




home and community care

A JOINT COMMONWEALTH AND STATE/TERRITORY PROGRAM
PROVIDING FUNDING AND ASSISTANCE FOR AUSTRALIANS IN NEED

Guidelines to assist you when

Choosing an Adjustable Bed

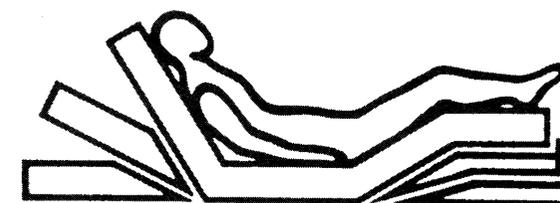

Adjustable Bed Options

Height Adjustment

This option is also referred to as a Hi / Lo option. It raises and lowers the mattress height to assist with:


- Getting in and out of bed easily
- Transferring to and from a wheelchair at the same height as the bed
- Making the bed
- Minimising the risk of injury to someone caring for a person in the bed
- Making transfers with a hoist easier

Head Elevation

This option raises the head of the bed, supporting the person's back and head. It can assist with:

- Sitting up in bed with support
- Getting out of bed
- Changing position in bed, to different angles

Leg Elevation

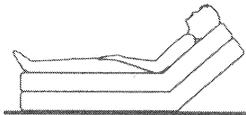

- Elevates the lower legs for comfort
- Assists with circulation and fluid drainage

Knee Break

- A leg elevation option
- Supports the thighs on an angle
- Keeps the lower legs elevated and parallel to the mattress
- Helps prevent a person from sliding down in the bed

Styles

- Ensemble-style beds with the above features but usually fixed in height
- Traditional hospital-style adjustable beds with wooden or metallic frames



The hospital style beds have space underneath for a hoist and are more suitable if 'monkey bars' are required. Many have the option of foot and head boards, giving them a more domestic look.

Methods of Operation

Manual or mechanical

- The height adjustment is often operated by an attendant with a winding handle or foot pedal
- The head elevation is manually adjusted by an attendant
- These tend to be the cheapest option but require more effort and assistance from others

Hydraulic

- Usually operated with a foot pump

Electric

- Operated by push-button or rocker switches for all adjustments
 - Modifications can be made for people who have difficulty pressing buttons
-